




## HANDSET INSTRUCTIONS IONIAN

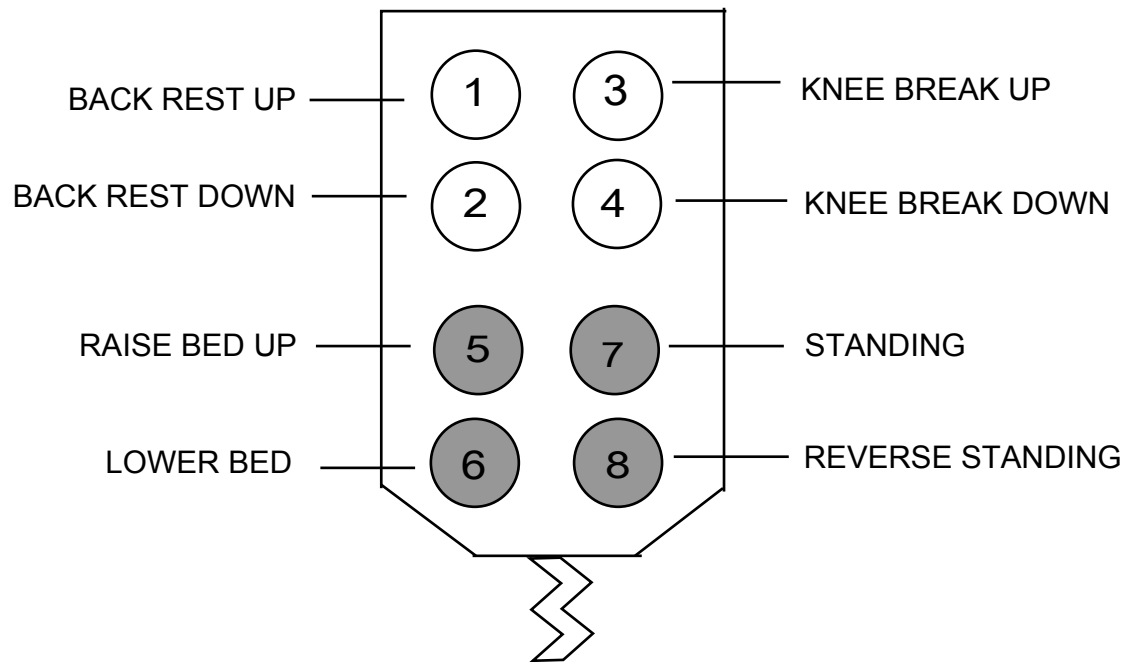

THE IONIAN STANDING BED SYSTEM IS OPERATED BY AN 8 BUTTON HANDSET.

1. BACK REST UP
2. BACK REST DOWN
3. KNEE BREAK UP
4. KNEE BREAK DOWN
5. STANDING BUTTON - THIS WILL ONLY OPERATE WHEN THE BED IS IN THE LOWERED POSITION. IF BED WON'T STAND PRESS BUTTON 6 UNTIL THE BED IS COMPLETELY LOWERED.
6. REVERSE STANDING - THIS BUTTON WILL LOWER THE BED TO FLAT POSITION.
7. RAISE BED UP
8. LOWER BED DOWN